

premio 30 evo

Laser duo

NEW!

Stimulate
the Qi in
depth



ACUPUNCTURE

AURICULOTHERAPY


ADDICTION

REFLEXOTHERAPY

Sedatelec[®]
YOUR WAY OF CARE,
OUR QUALITY DEVICES

premi^o 30 ev^o

Laser duo

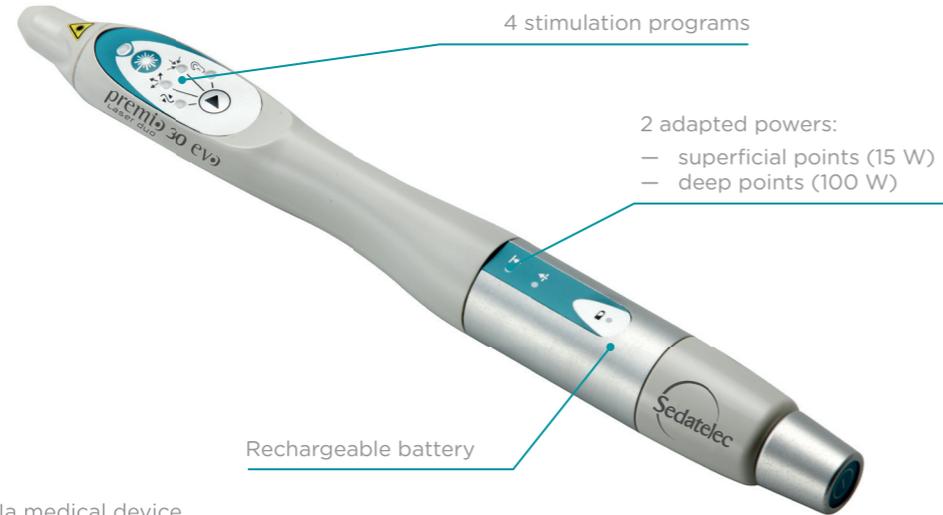

Class IIa medical device

A cold laser for painless infrared stimulation of Acupuncture, Auricular and Trigger points



Blending tradition and modernity: safe and effective

CORDLESS: easy to use, without restrictions

ECONOMICAL: no consumables

ERGONOMIC: direct access to 4 programmes via the dedicated button

ADAPTED: 2 adapted powers:

- Superficial points - 15 W: face, ears, children, etc.

- Deep points - 100 W: other body points

DYNAMIC: mobilises the Qi

FAST: 30 sec / point



Acupuncture stimulation

Harmonisation

This mode re-establishes meridian balance in the point and circulation by alternating the pulse rates.

Sedation

In just a few seconds, the Qi is dispersed and stagnant energy evacuated.

Tonification

Fast and controlled stimulation of the acupuncture points.

Auriculotherapy stimulation

Nogier frequency* sweep

Successively emitting all of the Nogier frequencies (from A to G), delivers an adapted message regardless of the symptomatology.

This sweep allows you to stimulate pathological auriculotherapy points, known for their effect on pain, psychological and functional disorders, addiction, and more...

Pulsed emissions for a dynamic message:

- Frequency wobulation to dynamise the Qi
- Nogier frequency* sweep in auriculotherapy

The Qi dynamic




The Qi expresses this idea of vital breath, energy and force in action, representing life and symbolising movement. This energy is continuously circulating in the meridians, expressing the balance of Yin and Yang and their regularly alternating cycle. It ensures that the living being functions harmoniously, as an individual integrated in both the immediate environment and the universe as a whole.

During the overall and in-depth clinical examination, after disorders have been determined and the causes analysed, **tonification** must be applied where there is a lack of energy, it must be **sedated** when there is an excess, and **harmonised**, with the flows reoriented when the crossing points are disturbed.


Harmonisation alternating tonification and sedation

Alternating between an acceleration and slowing of the emission rate.

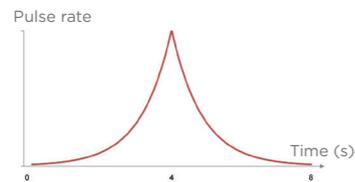


Sedation (10,000 to 100 Hz)

Obtained by slowing the emission pulse rate.

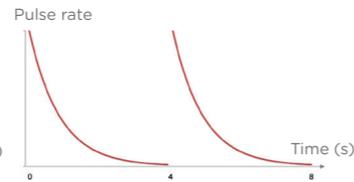

Tonification (100 to 10,000 Hz)

Acceleration of the pulse emission rate, for a dynamising effect.

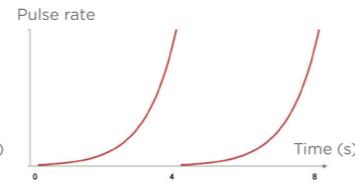

Wobbling to maintain the Qi:


Sedatelec has developed a specific approach to physiological impulses no longer based on repetitive frequencies but on dynamics, acceleration or deceleration as variation in the signal is considered to form the basis of the information decoded by the body.

What stimulates an organism is the dynamic, the change and the variation.

Effective stimulation


Wavelength

Infrared emission at a wavelength of 905 nm passes through the cutaneous barrier to penetrate deep cells safely, as it is less absorbed by water, haemoglobin and melanin.

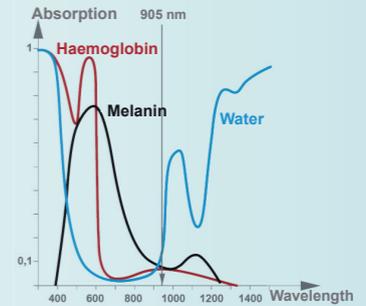

Peak power


You can set the power to **15** or **100 Watts** depending on the depth you wish to reach. Peak power of **100 Watts** gives a laser beam **with the power needed to penetrate deep cells** and impart the required information to the point.

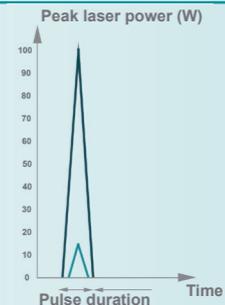

Pulsed emission


Very short (<100 ns) pulsed emission repeated with significant rest time. This imparts energy between:

- **150 mJ/30 sec** for 15 Watts
- **1.25 J/30 sec** for 100 Watts, with no destructive thermal effect.

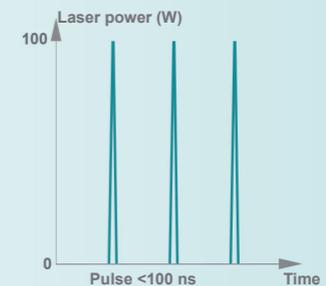



40 years of cold laser experience

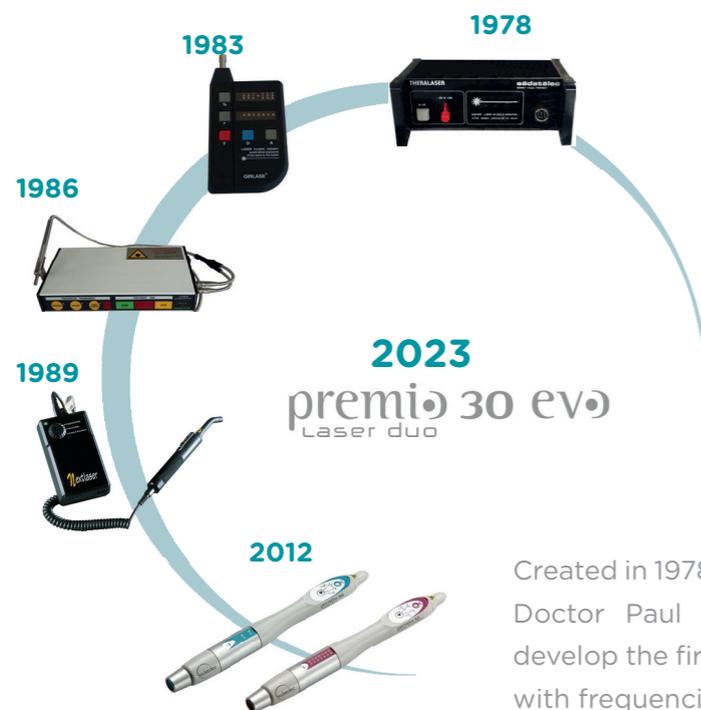

Created in 1978, Sedatelec worked with Doctor Paul Nogier (1908-1996), to develop the first laser-emitting devices with frequencies.

From R&D to marketing, everything, including assembly, is Made in France.

Sedatelec products are guaranteed for 2 years, with After-Sales Service provided at the plant in Lyon for up to 10 years after the product is removed from the market.

CONÇU & FABRIQUÉ
EN FRANCE



Soft Laser effects observed by all

– REFLEXOLOGY D.F. Lyon

The Premio 30 EVO is my reflexology, auricular and facial assistant. It is great for painless, easy to use stimulation. I use it to treat tendonitis and sciatic pain by stimulating the corresponding auricular points, and the effect felt is immediate.

I can no longer do without it!

– ACUPUNCTURE E.S. Aix-en-Provence

I've tested the Premio 30 EVO and I'm impressed by the results in my patients!

I use it to stimulate the acupuncture points using the "superficial" and "deep" modes. The effectiveness is genuinely comparable to that of a needle. I strongly recommend it!

– ACUPUNCTURE Nancy

My colleagues and I tested the PREMIO 30 EVO vs the last-generation PREMIO 30.

I myself am very enthusiastic, and my colleagues have also greatly appreciated it. We use it every day.

– TCM PRACTITIONER S.C. St Genis Laval

I have been using the PREMIO 30 for a year and am completely satisfied. It is a superb tool, specifically adapted to people who are fragile (children, seniors, etc.) and people with a fear of needles.

I recommend it.

testimony

A cordless Laser, ideal for both the surgery and house calls

The case includes:

- The PREMIO 30 EVO
- 1 emission control tester
- 1 pair of protective goggles
- 1 charger
- 1 user manual



Prior to any use, carefully read the user manual



*Laser device, IEC 60825-1 class 3B
Invisible laser radiation*

Exposure to the beam can harm the eyes.

Caution: Federal law restricts this device to sale by or on the order of a physician.

distributor's stamp

CE 0123
Made in France



Chemin des Mûriers • 69540 Irigny (Lyon) • France
Tel: +33 (0) 472 663 322 / Fax: +33 (0) 478 508 903
contact@sedatelec.com • www.sedatelec.com